
INGOLD

Leading Process Analytics

InPro 6950 Optical O₂ Sensors

Instruction manual Bedienungsanleitung Instructions d'utilisation

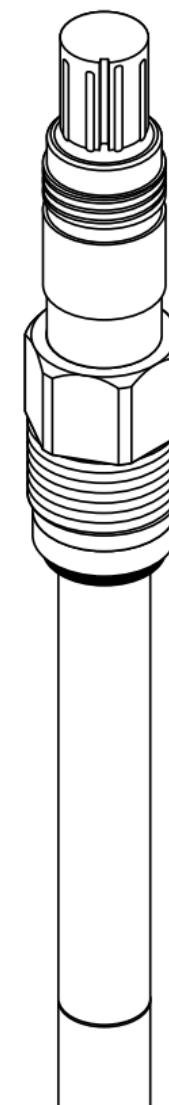InPro 6950 Optical O₂ Sensors

52 206 137 F

METTLER TOLEDO

English **Page** **3**

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InPro 6950 Series O₂ Sensors

Instruction manual

InPro, InFit, InTrac and ISM are trademarks of the METTLER TOLEDO Group.

Contents

1	Introduction.....	5
2	Important notes	6
2.1	Notes on operating instructions	6
2.2	Intended use	6
2.3	Safety instructions.....	7
2.4	Environmental protection	9
2.5	Examples of some typical applications.....	9
2.6	Use in Ex-zones	9
2.7	Ex-classification ATEX	10
2.7.1	Introduction.....	10
2.7.2	Rated data	10
2.7.3	Special conditions	11
2.8	Ex-classification – FM approved	12
2.8.1	Introduction.....	12
2.8.2	Rated data	13
2.8.3	Special conditions	13
2.8.4	Applied standards:	14
3	Product description.....	15
3.1	General information	15
3.2	Principle	15
3.3	Scope of delivery	16
3.4	Equipment features	17
4	Installation.....	18
4.1	Mounting the sensor	18
4.2	Connection	19
4.2.1	Connecting the InPro 6950 to a VP cable	19
4.2.2	Connecting the VP cable to the transmitter	19
5	Operation.....	20
5.1	Start-up and polarizing	20
5.2	Calibration	21
5.2.1	Purpose of calibration	21
5.2.2	What you have to know for calibration	21
5.2.3	Single point calibration	22
5.2.4	Dual point calibration.....	22
6	Maintenance	23
6.1	Inspection of the sensor.....	23
6.1.1	Visual inspection	23
6.1.2	Testing the METTLER TOLEDO O ₂ sensor master InPro 6950	24
6.1.3	Testing the sensor via a transmitter	25
6.2	Changing the electrolyte, the membrane body or the interior body	26
7	Storage.....	29
8	Product specification	30
8.1	Certificates	30
8.2	Specifications.....	31
9	Ordering information	32
9.1	Sensors	32
9.2	Accessories.....	32
9.3	Spare parts	32
9.4	Recommended transmitters.....	33
9.5	Recommended housings	33
10	Theory of the polarographic sensor.....	34
10.1	Introduction.....	34
10.2	Principle of the design of an oxygen electrode	35
10.3	Parameters determining current	36
10.4	Polarization voltage	37
10.5	Temperature	37
10.6	Dependence on flow	37
10.7	Oxygen partial pressure – oxygen concentration	38

1 Introduction

Thank you for buying the **InPro™ 6950 sensor from METTLER TOLEDO.**

The construction of the InPro 6950 sensors employs leading edge technology and complies with safety regulations currently in force. Notwithstanding this, improper use could lead to hazards for the user or a third-party, and/or adverse effects on the plant or other equipment.

Therefore, the operating instructions must be read and understood by the persons involved before work is started with the sensor.

The instruction manual must always be stored close at hand, in a place accessible to all people working with the InPro 6950.

If you have questions, which are not or insufficiently answered in this instruction manual, please contact your METTLER TOLEDO supplier. They will be glad to assist you.

2 Important notes

2.1 Notes on operating instructions

These operating instructions contain all the information needed for safe and proper use of the InPro 6950 sensor.

The operating instructions are intended for personnel entrusted with the operation and maintenance of the sensors. It is assumed that these persons are familiar with the equipment in which the sensor is installed.

Warning notices and symbols

This instruction manual identifies safety instructions and additional information by means of the following symbols:

This symbol draws attention to **safety instructions and warnings of potential danger** which, if neglected, could result in injury to persons and/or damage to property.

This symbol identifies **additional information and instructions** which, if neglected, could lead to defects, inefficient operation and possible loss of production.

2.2 Intended use

METTLER TOLEDO InPro 6950 sensors are intended solely for inline measurement of the oxygen partial pressure in liquids and gases, as described in this instruction manual.

Any use of these sensors which differs from or exceeds the scope of use described in this instruction manual will be regarded as inappropriate and incompatible with the intended purpose.

The manufacturer/supplier accepts no responsibility whatsoever for any damage resulting from such improper use. The risk is borne entirely by the user/operator.

Other prerequisites for appropriate use include:

- compliance with the instructions, notes and requirements set out in this instruction manual.
- acceptance of responsibility for regular inspection, maintenance and functional testing of all associated components, also including compliance with local operational and plant safety regulations.
- compliance with all information and warnings given in the documentation relating to the products used in conjunction with the sensor (housings, transmitters, etc.).
- observance of all safety regulations governing the equipment in which the sensor is installed.

- correct equipment operation in conformance with the prescribed environmental and operational conditions, and admissible installation positions.
- consultation with METTLER TOLEDO Process Analytics in the event of any uncertainties.

2.3 Safety instructions

- The plant operator must be fully aware of the potential risks and hazards attached to operation of the particular process or plant. The operator is responsible for correct training of the workforce, for signs and markings indicating sources of possible danger, and for the selection of appropriate, state-of-the-art instrumentation.
- It is essential that personnel involved in the commissioning, operation or maintenance of these sensors or of any of the associated equipment (e.g. housings, transmitters, etc.) be properly trained in the process itself, as well as in the use and handling of the associated equipment. This includes having read and understood this instruction manual.
- The safety of personnel as well as of the plant itself is ultimately the responsibility of the plant operator. This applies in particular in the case of plants operating in hazardous zones.
- The oxygen sensors and associated components have no effect on the process itself and cannot influence it in the sense of any form of control system.
- Maintenance and service intervals and schedules depend on the application conditions, composition of the sample media, plant equipment and significance of the safety control features of the measuring system. Processes vary considerably, so that schedules, where such are specified, can only be regarded as tentative and must in any case be individually established and verified by the plant operator.
- Where specific safeguards such as locks, labels, or redundant measuring systems are necessary, these must be provided by the plant operator.
- A defective sensor must neither be installed nor put into service.
- Only maintenance work described in this operating instruction may be performed on the sensors.
- When changing faulty components, use only original spare parts obtainable from your METTLER TOLEDO supplier (see spare parts list, "Section 9.3").

- No modifications to the sensors and the accessories are allowed. The manufacturer accepts no responsibility for damages caused by unauthorised modifications. The risk is borne entirely by the user.
- When using cables for sensor connection not supplied and recommended by the manufacturer the following standards need to be considered:
 - US National Electrical Code® (ANSI/NFPA 70 [NEC®]), where applicable.
 - Canadian Electrical (CE) Code® (CEC Part 1, CAN/CSA-C22.1), where applicable.
- For guidance on US installations, see ANSI/ISA-RP12.06.01, Installation of Intrinsically Safe Systems for Hazardous (Classified) Locations.
- Control room equipment connected to intrinsically safe associated apparatus in the US should not use or generate more than the specified Um of the associated apparatus.
- Care must be taken during installation to avoid impacts or friction that could create an ignition source.
- Tampering and replacement with non-factory components may adversely affect the safe use of the system.
- Insertion or withdrawal of removable electrical connectors or modules is to be accomplished only when the area is known to be free of flammable vapors.
- **WARNING !**
INTRINSICALLY SAFE APPARATUS CAN BE A SOURCE OF IGNITION IF INTERNAL SPACINGS ARE SHORTED OR CONNECTIONS OPENED.
- **WARNING !**
DO NOT OPEN WHEN AN EXPLOSIVE ATMOSPHERE IS PRESENT.
- **WARNING !**
SUBSTITUTION OF COMPONENTS MAY IMPAIR INTRINSIC SAFETY.
- **WARNING !**
SUBSTITUTION OF COMPONENTS MAY IMPAIR SUITABILITY OF THE EQUIPMENT.
- **WARNING !**
FOR CONNECTION ONLY TO NON-FLAMMABLE PROCESSES.
- **WARNING !**
TO MAINTAIN THE ENCLOSURE IP66 RATING, THE CONNECTOR MUST BE FULLY ENGAGED.

2.4 Environmental protection

Waste electrical products should not be disposed of with household waste. Please recycle where facilities exist. Check with your Local Authority or retailer for recycling advice.

2.5 Examples of some typical applications

Below is a list of examples of typical fields of application for the oxygen sensors. This list is not exhaustive.

Measurement in liquids:

- Brewing
- Beverage filtration
- Filling stations
- Spring water conditioning

Measurement in gases:

- CO₂ recovery
- CO₂ purity
- Product storage
- Inert production

2.6 Use in Ex-zones

Attention!

For an installation in Ex-zones please read the guidelines following hereafter:

Ex-classification ATEX:

- Ex II 1/2G Ex ia IIC T6...T3 Ga/Gb
Ex II 1/2D Ex ia IIIC T69 °C...T161 °C Ga/Db

Number of the test certificate:

SEV 14 ATEX 0169X
IECEx SEV 14.0026X
CML 22 UKEX 2212X

Ex-classification FM approved:

IS / I, II, III / 1 / ABCDEFG / T6 Ta = 60 °C
- 53 800 002; Entity

Project ID: 3021227

FM certificate number: FM16US0034X

FM18CA0021X

2.7 Ex-classification ATEX

2.7.1 Introduction

According to Directive 2014/34/EU (ATEX 114)¹⁾ Appendix I, InPro 6XXX*/**/*/* oxygen sensors are devices group II, category 1/2G and according to RL 99/92/EG (ATEX 137)²⁾ may be used in zones 0/1 or 1/2 or 1 or 2 and gas groups IIA, IIB and IIC that are potentially explosive due to combustible substances in above stated temperature classes.

For use / installation, the requirements of EN 60079-14 must be observed.

According to Directive 2014/34/EU (ATEX 114)¹⁾ Appendix I, InPro 6XXX*/**/*/* oxygen sensors are devices group II, category 1/2D and according to RL 99/92/EG (ATEX 137)²⁾ may also be used in zones 20/21 or 21 or 22 that contain combustible dusts.

The digital sensor circuit is part of a common intrinsically safe system and is for operation connected to a separately certified transmitter.

The digital sensor circuit as part of an intrinsically safe system is isolated from the not-intrinsically safe electric circuits up to a maximum rated voltage of 375 V and from grounded parts up to a maximum rated voltage of 30 V.

2.7.2 Rated data

Measuring circuit

Intrinsically safe ignition protection class Ex ia IIC only for connection to a certified intrinsically safe circuit.

Analog O₂ oxygen sensor maximum values

$U_i \leq 16$ V

$I_i \leq 190$ mA

$P_i \leq 200$ mW

$L_i = 0$ (effective internal inductance)

$C_i = 900$ pF (effective internal capacitance)

 **Note:** The above maximum values are each the total of all individual circuits of the associated intrinsically safe power supply and transmitter.

¹⁾ For UK statutory Requirements SI 2016 No. 1107

²⁾ For UK regulations 7 and 11 of the Dangerous Substances and Explosive Atmospheres Regulations 2002 (DSEAR)

2.7.3 Special conditions

- The relationship between the maximum permissible ambient or media temperature and temperature class, for category 1G applications, zone 0, is shown in the following table:

Temperature class	Max. ambient or media temperature
T 6	68 °C
T 5	80 °C
T 4	108 °C
T 3	160 °C

- The relationship between the maximum permissible ambient or media temperature and temperature class, for category 1D applications, zone 20, is shown in the following table:

Temperature class	Max. ambient or media temperature
T 69 °C	68 °C
T 81 °C	80 °C
T 109 °C	108 °C
T 161 °C	160 °C


- The capacitance and inductance of the connecting cable has to be considered.
- The O₂ Oxygen sensor type InPro 6XXX can be used in / with the fittings InFit™ 76*-*** or InTrac™ 7**-***, or in / with other suitable fittings in potentially explosive areas.
- The metal body of the O₂ Oxygen sensors, or the fittings InFit 76*-*** or InTrac 7**-***, or other appropriate fitting is optionally included in the routine pressure test of the system.
- The independent fitting used for installation of O₂ Oxygen sensor must be conductively connected to the equipotential bonding system.

2.8 Ex-classification – FM approved

2.8.1 Introduction

The following FM control drawing and the standards listed in section 2.8.4 must be observed, where applicable:

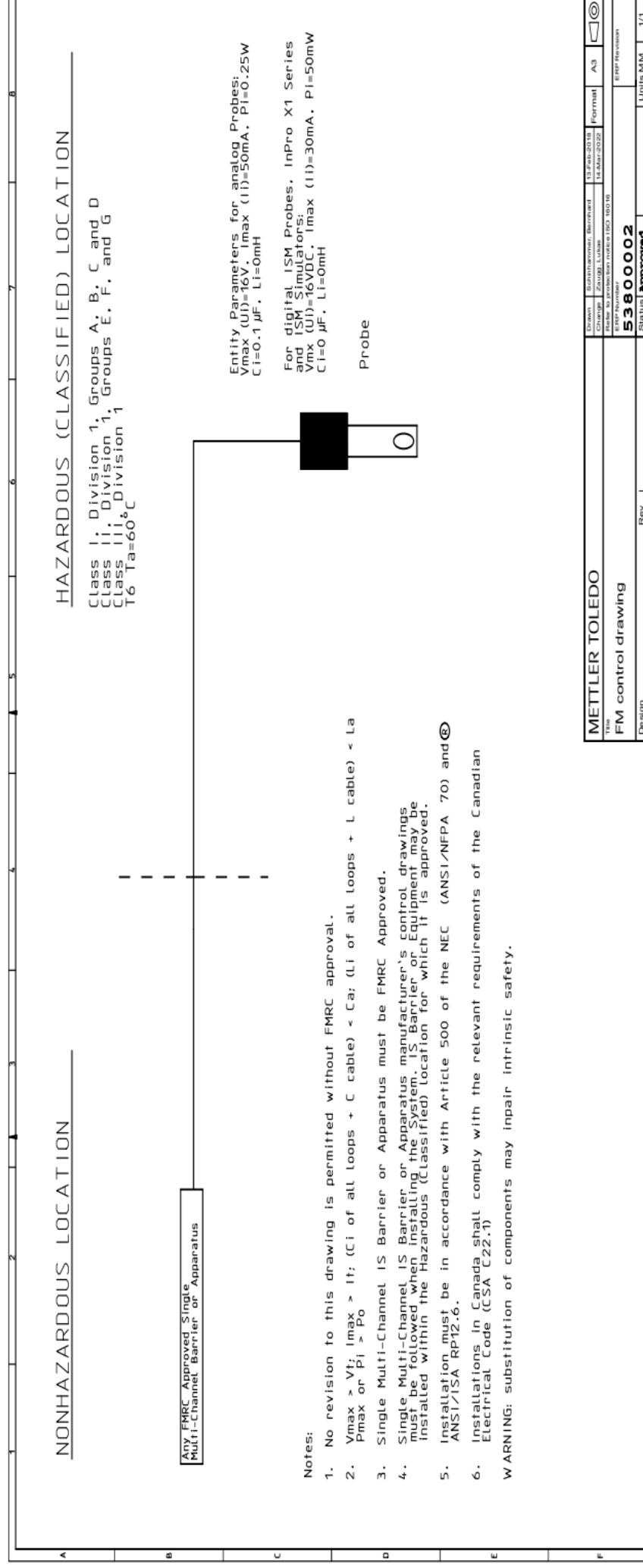2.8.2 Rated data

For rated data please observe section 2.7.2.

2.8.3 Special conditions

In type of protection intrinsically safe apparatus, the analytical probes (InPro 6000 Series with InFit 76X/Y and InTrac 7XX/YY Series housings) equipment is designated with the following specific conditions of use.

1. The dissolved oxygen sensor shall be installed in compliance with the enclosure, mounting, spacing and segregation requirements of the ultimate application, including a tool removable cover, and is suitable for use with FM Approved InFit 76X/Y and InTrac 7XX/YY Series housings.
2. Process temperature no greater than +130 °C.
3. Maximum permissible working pressure is 12 barg (174 psig).
4. Potential Electrostatic Charging Hazard – To prevent the risk of electrostatic sparking, the non-metallic surface should only be cleaned with a damp cloth.
5. Enclosures containing titanium constitute a potential risk of ignition by impact or friction. Care must be taken into account during installation and use to prevent impact or friction.

2.8.4 Applied standards:

United States Standards

Title	Number	Issue Date
Approval Standard for Electrical Equipment for Use in Hazardous (Classified) Locations – General Requirements	FM Class 3600	2022
Approval Standard for Intrinsically Safe Apparatus and Associated Apparatus for Use in Class I, II & III, Division 1, Hazardous (Classified) Locations	FM Class 3610	2021
Approval Standard for Electrical Equipment for Measurement, Control and Laboratory Use	FM Class 3810	2021
Explosive Atmospheres – Part 0: Equipment – General Requirements	ANSI/UL 60079-0:2020	2020
Explosive Atmospheres – Part 11: Equipment Protection by Intrinsic Safety "i"	ANSI/UL 60079-11:2018	2018
Safety Requirements for Electrical Equipment for Measurement, Control, and Laboratory Use – Part 1: General Requirements	ANSI/UL 61010-1:2018	2018

Canadian Standards

Title	Number	Issue Date
Explosive Atmospheres – Part 0: Equipment – General Requirements	CAN/CSA-C22.2 No. 60079-0:2019	2019
Explosive Atmospheres – Part 11: Equipment Protection by Intrinsic Safety "i"	CAN/CSA-C22.2 No. 60079-11-14:R2018	R2018
Safety Requirements for Electrical Equipment for Measurement, Control, and Laboratory Use – Part 1: General Requirements	CAN/CSA-C22.2 No. 61010-1-12:R2017	R2017

3 Product description

3.1 General information

The **oxygen sensor InPro 6950** with integrated temperature probe is used for measurement of oxygen at low concentrations.

The sensor is **sterilizable** and **compatible with CIP** (cleaning in place).


InPro 6950 sensors with ISM™ functionality offer Plug and Measure as well as enhanced diagnostics features.

3.2 Principle

The InPro 6950 is a Clark polarographic O₂ sensor that has been expanded with a separate reference and guard ring electrode. This 4-electrode system is characterized as follows:

- a) The measurement system of this InPro 6950 sensor consists of a working electrode (cathode), a counter electrode (anode), a reference electrode and a guard ring electrode. The measurement system is separated from the process medium by an oxygen permeable membrane.
- b) The transmitter supplies a constant voltage to the cathode and the guard ring electrode so that oxygen can be reduced on both electrodes.
- c) The oxygen molecules migrate from the measurement medium through the membrane to the cathode to which the voltage is supplied and are reduced. At the same time, oxidation takes place at the anode and oxygen and water are released into the electrolyte. This causes a current to flow between the anode and cathode which is directly proportional to the partial pressure of oxygen (pO₂) in the process medium.
- d) The guard ring electrode built into the sensor reduces the oxygen that migrates to the side of the cathode and can distort the measurement. The guard ring electrode therefore enables traces of oxygen to be precisely determined at even the lowest concentrations.

The 4-electrode system of the InPro 6950 guarantees a high degree of precision, rapid response and a low detection limit.

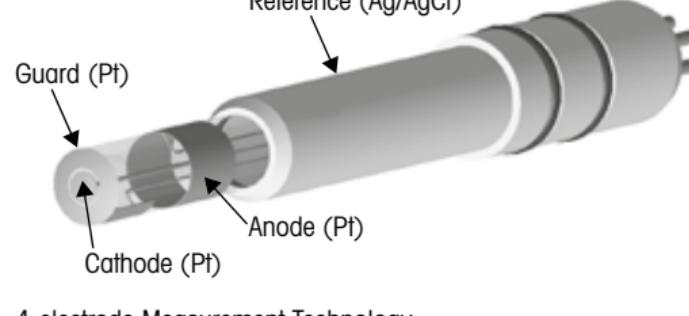4-electrode Measurement Technology

Please refer to "Section 10 – Theory of the polarographic sensor" for further information.

3.3 Scope of delivery

Each InPro 6950 sensor is supplied fully assembled and factory-tested for correct function together with:


- a quality control certificate
- inspection certificates 3.1
(complying with EN 10204)

METTLER TOLEDO DO sensors are supplied with fitted membrane body filled with electrolyte and covered with the watering cap and have been checked for proper function.

See electrolyte order details at section „Spare Parts”, 9.3.

3.4 Equipment features

12 mm sensor

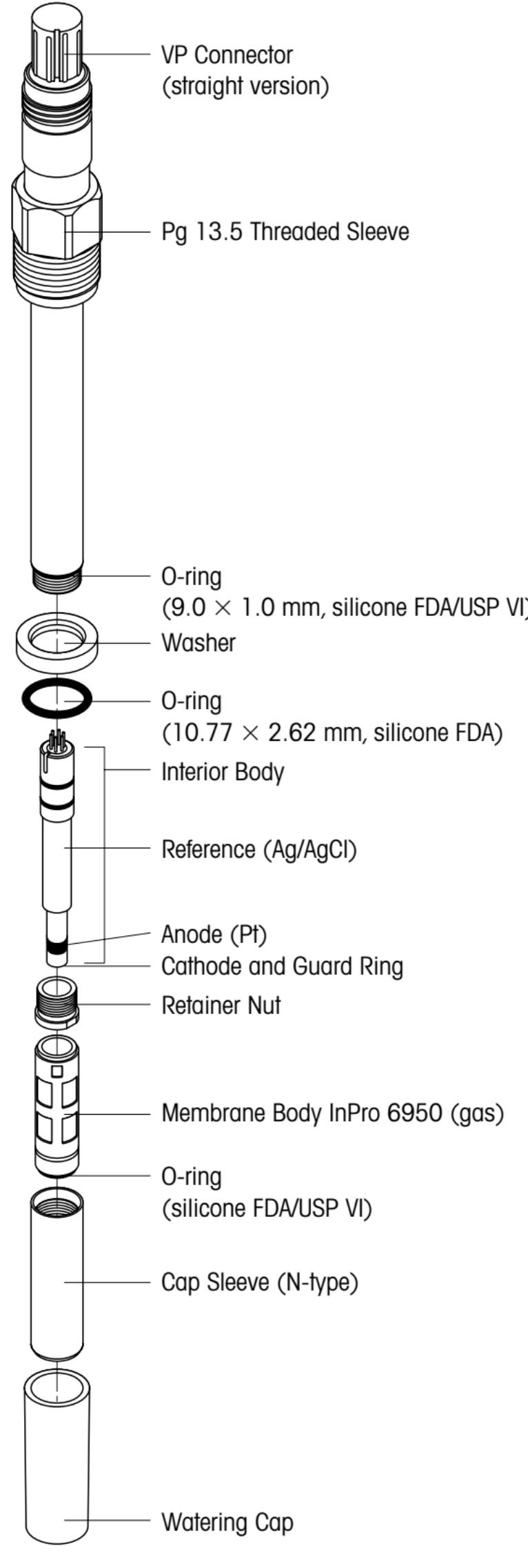4 Installation

4.1 Mounting the sensor


Important! Remove the watering cap before mounting the sensor.

Mounting the sensor in a housing

Please refer to the instruction manual of your housing explaining on how to mount the sensor in place.

Mounting the sensor directly on a pipe or a vessel

The **12 mm sensors** can be mounted directly through a socket with inside thread Pg 13.5 and securely tightened via the Pg 13.5 threaded sleeve.

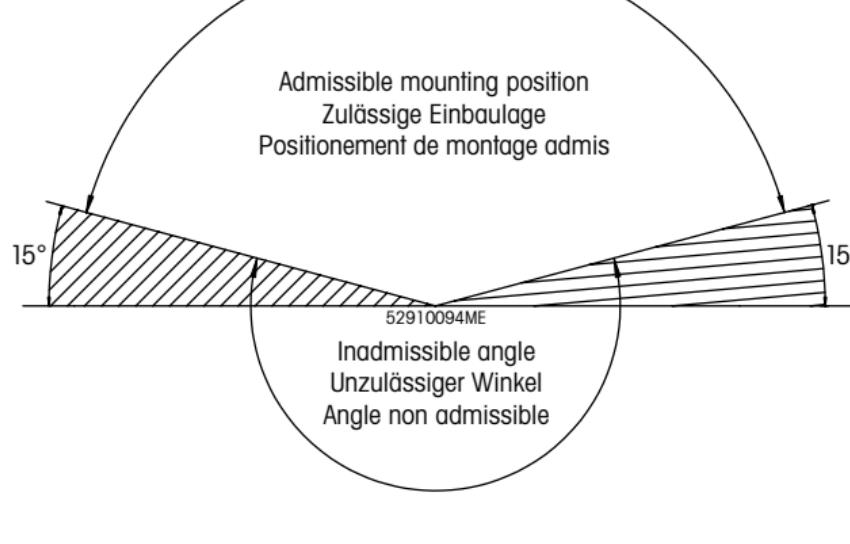
4.2 Connection

4.2.1 Connecting the InPro 6950 to a VP cable


 The sensor is connected to the transmitter via a VP cable. The VP cable ensures a secure connection between the transmitter and the sensor under harsh industrial conditions. The robust watertight IP 68 connector housing guarantees maximum process safety.

 **Do not touch the sensor at the VP connector plug!**

To connect the VP cable to the sensor align the slit of the VP connector with the pin in the plug. Then tightly screw the plug to fasten the two parts.

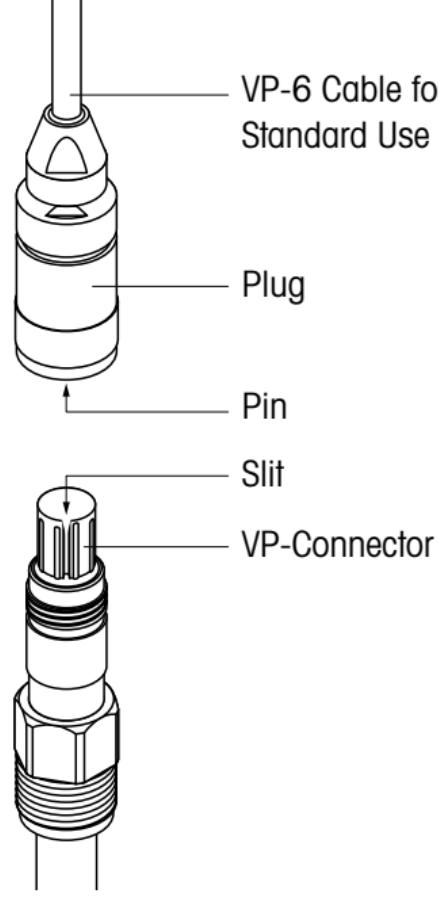4.2.2 Connecting the VP cable to the transmitter

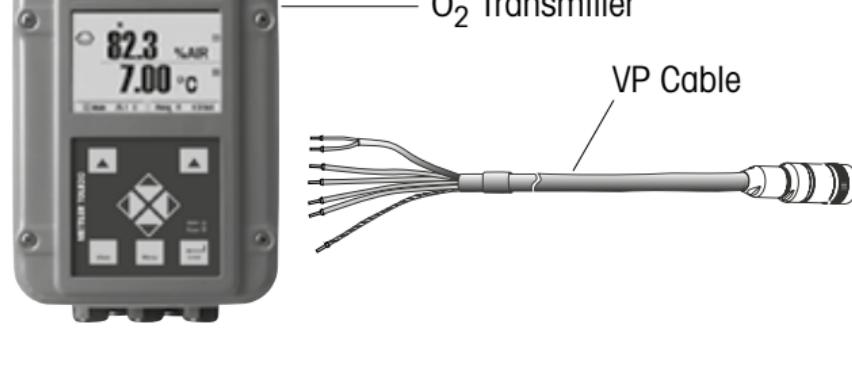 **Note:** Cable assignment can be found in the METTLER TOLEDO VP cable instruction manual.

 **Note:** For connecting the cable to the terminals of the transmitter, please refer to the instructions given in the METTLER TOLEDO transmitter manual.

5 Operation

Important! Before using the sensors for the first time, the electrolyte should be replaced (see "Chapter 6.2"). Due to possible adverse conditions during transport and storage (e.g. airfreight; pressure and temperature variations), the quality of the electrolyte may become impaired. Poor electrolyte quality can lead to erroneous measurement values.

5.1 Start-up and polarizing

Important! The watering cap must be removed before mounting the sensor in the process.

When the system is operated for the first time or if the sensor has been disconnected from the voltage source (transmitter or O₂ sensor master) for longer than 5 minutes, the sensor has to be polarized prior to calibration by connecting it to the operating O₂ transmitter or to a sensor master. After 6 hours, the sensor is fully polarized and ready for operation. During this time, it is recommended to keep the watering cap filled with conditioning solution on the sensor (see "Section 7 – Storage"), especially if the sensor is being polarized for longer than 6 hours.

A shorter polarization period is sufficient if the sensor has been disconnected for only a few minutes. The following table serves to establish the correct polarization time in relation to the depolarization time.

Depolarization time ¹ t_{depol} [Min.]	Minimum required polarization time ² [Min.]
$t_{\text{depol}} > 30$	360
$30 > t_{\text{depol}} > 15$	$6 \times t_{\text{depol}}$
$15 > t_{\text{depol}} > 5$	$4 \times t_{\text{depol}}$
$t_{\text{depol}} < 5$	$2 \times t_{\text{depol}}$

1 Depolarization time: Time span in which the polarization voltage is cut off from the sensor. This is the case during:

- change of electrolyte
- change of membrane body
- the time the cable is disconnected or no transmitter or sensor master is connected to the cable.

2 Polarization time: Time span during which the sensor is under a polarization voltage.

Important! Setting of the polarization voltage on the transmitter for correct measurements:

- Standard applications: – **675 mV**

- Measurement of permanently low oxygen concentrations (< 500 ppb in liquids or < 10,000 ppm [vol.] in gases) in the presence of volatile acidic components (e.g. carbon dioxide during measurements in breweries): – **500 mV**

 **Note:** To ensure the supply of the correct polarization voltage the transmitter must be set accordingly.

5.2 Calibration

5.2.1 Purpose of calibration

Each oxygen sensor has its own individual slope and own individual zero point. Both values are subject to change, for example, through electrolyte consumption or after exchange of electrolyte or membrane body. To ensure high measurement accuracy of the sensor, a calibration must be carried out after each change of electrolyte or membrane. Prior to calibration, the sensor has to be polarized for at least 6 hours.

 **Please remove the watering cap from the sensor, rinse the sensor with water and dry it for at least 10 minutes before starting calibration.**

 **To check if your sensor needs a recalibration, you may dry it and take it in the air to check that the reading is close to 100 %. If not, then the sensor needs a new calibration.**

5.2.2 What you have to know for calibration

Calibration is generally carried out at the polarization voltage at which measurements are taken.

Zero-point calibration is recommended:

- when measuring less than 5 ppb in liquids
- when measuring less than 125 ppm (by vol.) in gases

The automatic calibration routine in the M 700 transmitter (SW700-011 "CO₂ compensation") is not suitable for the InPro 6950 sensor.

 **General remarks:**

- **For calibration in air, the sensor membrane must be dry**, since adhering water drops can falsify the measured oxygen value.
- Make sure that the **oxygen saturation index** of the calibration medium is **correct and remains constant** during calibration.
- In the event of calibration in water or sample medium, the **calibration medium must be in equilibrium with the air**. Oxygen exchange between water and air is only very slow. Therefore it takes quite long time until water is saturated with atmospheric oxygen.

- For correct calibration, a minimum flow rate of the calibration medium is necessary.
- Make sure that all other parameters, such as temperature and pressure, are constant.

For continuous applications, we recommend **periodic recalibration in line with your requirements on accuracy, the type of process in operation and your own experience**. The frequency of the need for re-calibration depends very much on the specific application, and therefore appropriate intervals cannot be exactly defined here.

5.2.3 Single point calibration

By carrying out a single point calibration, the factual slope of the sensor can be established. The calibration medium can be water with known oxygen saturation index (e.g. air-saturated water) or air with known water-vapor saturation (e.g. water-vapor saturated air).

After the sensor signal has stabilized, the complete measuring system can then be calibrated to the 100 % value of the desired measurable variable, e.g. 100 % air, 20.95 % O₂, or 8.26 ppm at 25 °C (77 °F) and normal pressure (see instruction manual for the transmitter).

5.2.4 Dual point calibration

By carrying out a dual point calibration both slope and zero point of the sensor can be established.

Important! In case of a **dual point calibration**, **always start by the zero point calibration** before calibrating the slope.

As a rule, the zero point should be adjusted to zero manually, or it is automatically performed by the transmitter (see instruction manual for the transmitter).

Zero-point calibration is recommended:

- when measuring less than 5 ppb in liquids
- when measuring less than 125 ppm (by vol.) in gases

Attention! Incorrect zero point calibration is a frequent source of measurement error. For correct calibration, we recommend the use of carbon dioxide gas with a level of purity of at least 99.9995 %.

After the sensor signal has stabilized (after 6 ... 12 hours), the sensor can be calibrated through the relevant transmitter to the 100 % value of the desired measurable variable, e.g. 0 % air, 0.0 % O₂, or 0.0 ppm (see instruction manual for the transmitter).

6 Maintenance

6.1 Inspection of the sensor

6.1.1 Visual inspection

To check your sensor, we recommend the following procedure:

- The contacts of the connector must be dry. Moisture, corrosion and dirt in the connector can lead to false readings.
- Check the cable for buckling, brittle areas or ruptures.
- Before calibration always examine the membrane foil optically for signs of damage. The foil must be intact and clean. Dirty membranes should be wiped clean using a soft, moist tissue.

Note: An undulated membrane has no influence on the sensor performance, assuming the membrane is intact.

- The membrane body must be replaced if the sensor has too long a response time, the reading is unstable or subject to drift, and if the sensor cannot be calibrated or the membrane shows sign of mechanical damage.
- Check the cathode area for discoloration, contamination or cracks in the glass. If necessary rinse with demineralized water and clean with a clean soft brush or soft paper tissue.

Attention! Do not use any cleaning agents containing alcohol. This could damage the sensor or lead to fault current.


Attention! The glass body is fragile and sensitive to vibration.

6.1.2 Testing the METTLER TOLEDO O₂ sensor master InPro 6950

We recommend to use the METTLER TOLEDO O₂ sensor master InPro 6950 to check the quality of your sensor as follows:

- Connect the sensor to the O₂ sensor master.

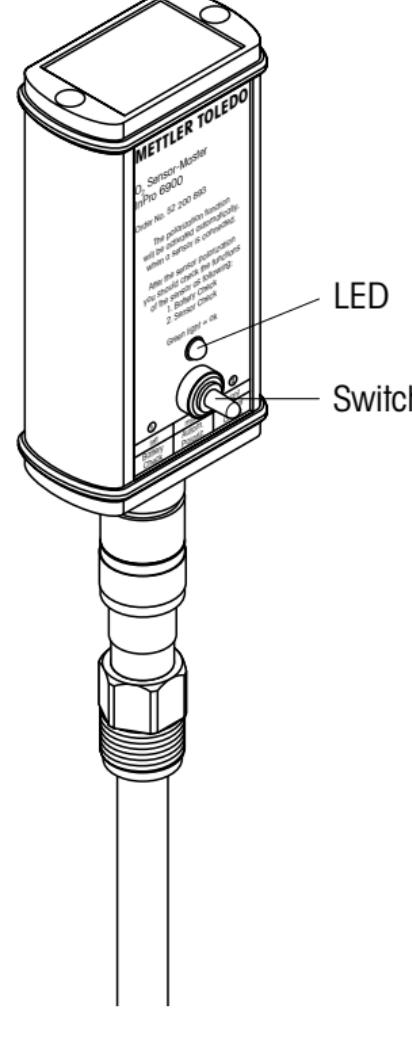As soon as the sensor is connected to the O₂ sensor master, the polarization function is automatically activated. Please note: if the sensor was disconnected from the transmitter for longer than 5 minutes, the sensor must be polarized first (polarizing time see "Section 5.1") to get representative test results.

- **Battery Check:**

Push the switch to the left. If the battery is ok and the O₂ sensor master is operational the green LED lights up. Otherwise, please consult the instruction manual of the O₂ sensor master.

- **Sensor Check:**

For this test the O₂ sensor must be **fully polarized** and the **membrane of the sensor must be dry and clean**.

Expose the sensor (connected to O₂ sensor master) to the air. By pushing the switch to the right to the position "2 – Sensor check", the O₂ sensor master checks whether the electrode current for air measurement delivered by the sensor is within the admissible range (2500 to 6000 nA for InPro 6950).

If the green LED lights up the current for air measurement is within the admissible range.

If the LED does not light up, you should check the battery of the O₂ sensor master (see instruction manual "Accessories"). If the battery is working, than there is probably a problem with your sensor. You should change the electrolyte and/or the membrane body of your sensor. If after a membrane change the LED still does not light up, this means that there is maybe something wrong with the interior body of the sensor. You should then change it (see "Section 6.2").

Important! The Sensor Check function only verifies the correctness of the electrode current for air measurement. In order to be absolutely sure of the functionality of the sensor, the residual signal in an oxygen free medium should also be controlled (see "Section 6.1.3").

Important! The InPro 6950 sensor is not compatible with the O₂ sensor master of the InPro 6900 or any other conventional polarization unit.

The 4-electrode system of the InPro 6950 sensor can only be checked and/or polarized directly on the transmitter or with the InPro 6950 O₂ sensor master.

6.1.3 Testing the sensor via a transmitter

A periodic zero current measurement (**no zero point calibration!**) is recommended for verification of proper sensor function.

Note: At the time you carry out the zero current measurement, the sensor must be polarized.

Zero current measurement can be done by using zeroing gel (order no. 30 300 435) or nitrogen (N₂) or carbon dioxide (CO₂) calibration gases with a purity of at least 99.995 %, alternatively in a sample medium saturated with one of these gases.

After 2 minutes in an oxygen-free sample medium, the reading on the transmitter should drop to below 10 % of the reading in ambient air, and within 10 minutes the value should have dropped to below 1 %.

If the measured values are too high, this suggests a depleted electrolyte or a defective membrane. In the first instance replace the electrolyte, and in the second case exchange both the membrane body and the electrolyte accordingly.

If after such procedures the above mentioned values are still not reached, replace the interior body. If this doesn't solve the problem too send the sensor to your local METTLER TOLEDO representative for inspection.

Many sample media contain volatile substances which, even at very low concentrations, have a clearly perceptible smell. Similarly to oxygen, these substances are able to invade the electrolyte through the gas-permeable membrane. Accordingly, they become noticeable when changing the electrolyte. In most cases, such substances have absolutely no influence on the measuring properties of the sensor. Slight discoloration of the electrolyte also has no effect on the measuring properties.

6.2 Changing the electrolyte, the membrane body or the interior body

 **Note:** The InPro 6950 uses a special electrolyte which contains an oxygen scavenger. This electrolyte ensures fast response time and together with the guard ring increases the signal stability of the sensor. The electrolyte must be exchanged on a regular basis or when the sensor has been exposed to the air for more than 24 hours without a watering cap or if the sensor has been stored for several months.

 **Note:** The content of the electrolyte bottle has to be used within 3 months after opening.

If the membrane and/or the interior body exhibits signs of failure (long response time, increased current in an oxygen-free medium, mechanical damage, etc.) the membrane body and/or the interior body has to be replaced.

 **Warning!** The O₂ electrolyte has an alkaline pH value.

Contact of electrolyte with mucous membrane or eyes is to be avoided. **Therefore protective gloves and safety glasses have to be worn for the following dismantling works.**

If such contact occurs, the affected area should be well rinsed with water. In the case of accident, or should ever any adverse signs appear, get immediate medical attention.

When changing the electrolyte, the membrane body or the interior body, please observe the following instructions (see also the following illustration):

Attention! Make sure that this maintenance step is carried out in **clean place**.

1. Unscrew the cap sleeve from the sensor shaft and carefully pull it off the sensor.
2. Pull off the membrane body from the interior body. If it is tight-fitted, eject by pushing it with the flat finger tip. Before electrolyte is refilled, the membrane body must be removed from the cap sleeve!

3. Rinse the interior body with demineralized water and carefully dab it dry with a paper tissue.


Note: steps 4 to 7 may only be carried out when changing the interior body.

4. Unscrew the retainer nut of the interior body with an adjustable wrench or with a 3/8" wrench.
5. Remove the interior body by pulling it out of the sensor shaft. If necessary use a plier.
Warning! Do not twist the interior body. Otherwise the connection pins can be damaged.
6. Insert the new interior body in the sensor shaft. Turn the interior body in the shaft until the slit of the interior body is aligned with the pin placed in the shaft.
7. Press the body in the shaft and screw the new retainer nut in place.
8. Examine the O-rings visually for mechanical defects, and replace if necessary.
9. Half-fill the new membrane body with O₂ electrolyte.

Note: The electrolyte bottle is equipped with a special pouring system. To ensure proper functioning, hold the bottle vertically, upside-down.

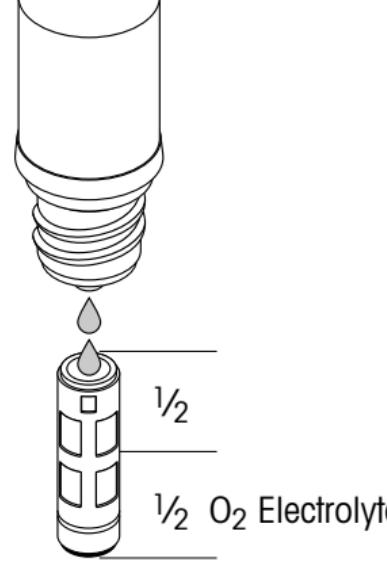Note: make sure that all air bubbles are removed from the membrane body. Air bubbles can be removed by carefully tapping on the membrane body.


10. Slip the membrane body over the interior body while holding the sensor in a vertical position. The excess electrolyte will be displaced and have to be removed with a paper tissue.

Important! No electrolyte, sample media or contamination may be present between the membrane body and the cap sleeve. Please check carefully!

11. Carefully slip the cap sleeve over the fitted membrane body, holding the sensor in a vertical position and screw it tight. The cap sleeve must be clean and dry.
12. After each exchange of electrolyte or membrane body, the sensor has to be repolarized and recalibrated.

Replacement of the interior body

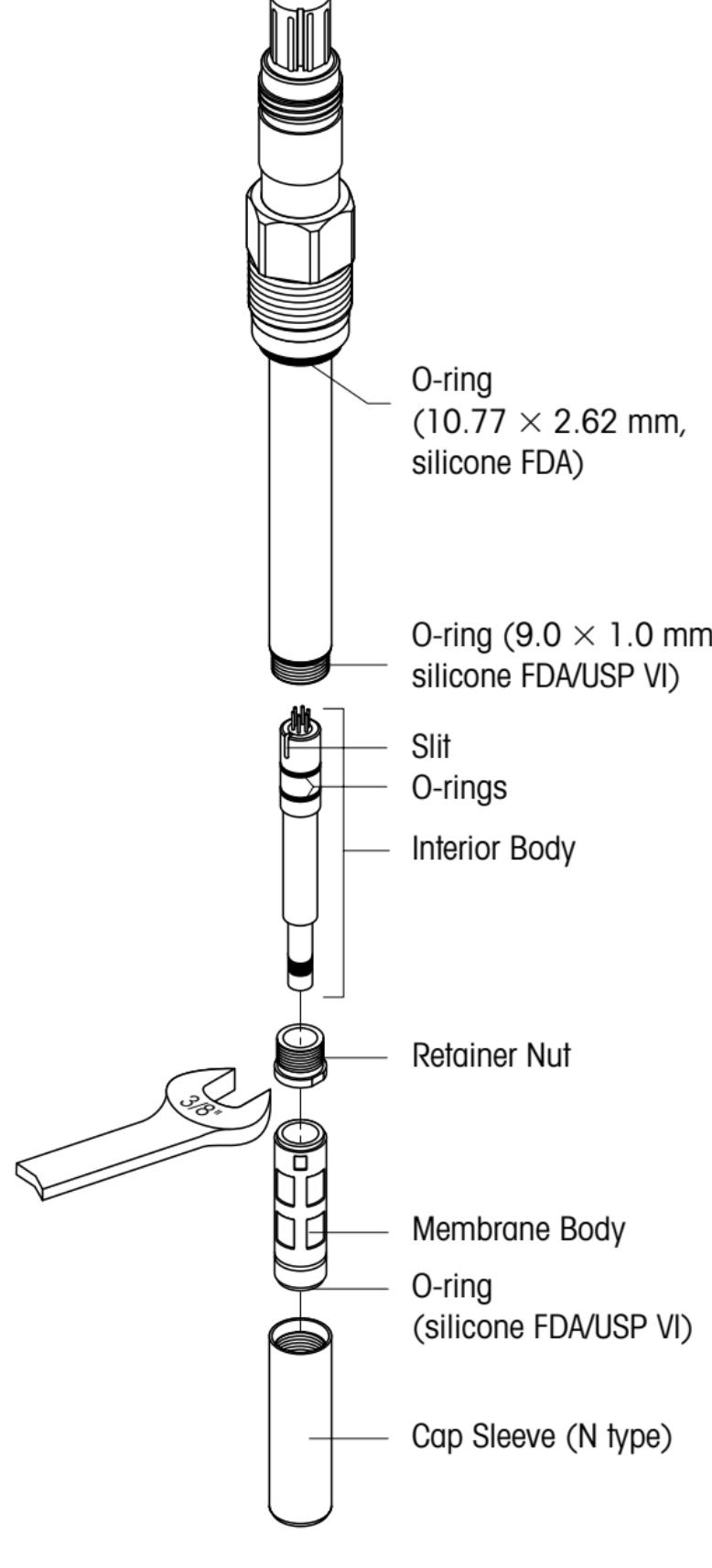7 Storage

For storage periods longer than 24 hours, we recommend you to use a watering cap filled with the same cleaning and conditioning solution as that used for our portable DO measurement system (order no. 52 200 255). This solution contains an oxygen scavenger which ensures that the electrolyte of the sensor does not come into contact with oxygen when not in use.

Prepare the cleaning and conditioning solution as follows: Dissolve one tablet in 40 ml of deionized water and wait 5 minutes for the tablet to be completely dissolved. Fill the watering cap with this solution and place it over the tip of the sensor. This solution has some cleaning properties which will keep the membrane free of microorganisms. In case you do not have any cleaning and conditioning set, you may also use checking gel or deaerated water in the watering cap. Before mounting the sensor in the process, always remove the watering cap and rinse the tip of the sensor with water.

Attention! The electrolyte should be replaced before use if the storage period exceeds 3 months.

Attention! If the storage period exceeds 6 months, the sensor should be **stored dry**, i.e. without any electrolyte in the membrane body. A sensor being stored dry (without electrolyte in the membrane body) **may on no account be connected to the O₂ sensor master InPro 6950 or any other polarization module.**

8 Product specification

8.1 Certificates

Each sensor is delivered with a set of **3.1 certificates** (complying with EN 10204).

All wetted metal parts (sensor shaft, cap sleeve and membrane body) are identified with an engraved symbol corresponding to the heat number on the paper certificate delivered with the sensor.

Each wetted metal part (sensor shaft, cap sleeve and membrane body) is polished in order to get a surface roughness lower than 0.4 µm (16 µin). This represents a roughness grade number of N5 (according to ISO 1320:1992).

8.2 Specifications

InPro 6950	
Measurement principle	Amperometric/Polarographic
Working conditions	
Pressure resistance measurement	0.2 ... 6 bar absolute [2.9 ... 87.0 psi absolute] gas applications
	0.2 ... 9 bar absolute [2.9 ... 130 psi absolute] liquid applications
Mechanical pressure resistance	Max. 12 bar absolute [174.0 psi absolute]
Measuring temperature range	0 ... 70 °C [32 ... 158 °F] gas applications
	0 ... 80 °C [32 ... 176 °F] liquid applications
Temperature range (mechanical)	– 5 ... 121 °C [23 ... 250 °F] (sterilizable)
Construction	
Temperature compensation	Automatic with built-in RTD
Cable connection	VarioPin (IP 68), straight or angled
O-ring material	Silicone FDA and USP VI approved
Membrane material	PTFE/Silicone
Wetted metal parts (3.1)	Stainless steel Special material on request
Surface roughness of wetted metal parts (ISO 1320:1992)	N5 ($R_a < 0.4 \mu\text{m}$ [$16 \mu\text{in}$])
Quick disconnect interior body	Traces
Cathode	Pt
Anode	Pt
Guard ring	Yes (Pt)
Reference	Ag/AgCl
Dimensions	
Sensor diameter	12 mm [0.47"]
Immersion length (a) for 12 mm sensor	70, 120, 220, 320 mm [2.8, 4.7, 8.66, 12.6"]
Performances	
Detection limit	0.1 ppb (water) 0.25 ppb (water / CO ₂) 5 ppm (vol.) gas
Accuracy	$\pm [1 \% + \pm 0.1 \text{ ppb}]$ of reading in liquids $\pm [1 \% + \pm 0.25 \text{ ppb}]$ of reading in CO ₂ containing liquids $\pm [1 \% + \pm 5 \text{ ppm}]$ of reading in gas
Response time at 25 °C/77 °F (air → N ₂)	98 % of final value < 90s
Sensor signal in ambient air (25 °C/77 °F)	2500 ... 6000 nA
Residual signal in oxygen-free medium	< 0.025 % of the signal in ambient air
Certification	
EHEDG, 3A	Yes (pending)
3.1 (EN 10204)	Yes
ATEX, IECEx and UKEX Certificate	Yes
FM Approval	Yes
FDA / USP VI	Yes
Quality Control	Yes
Compatibility	
with METTLER TOLEDO transmitters	see "Section 9.4"
with METTLER TOLEDO housings	see "Section 9.5"

9 Ordering information

For more detailed information refer to the technical data sheet. Ask your local distributor.

9.1 Sensors

Sensor	Order No.
InPro 6950/12/120/NP (4.7")	52 206 098

9.2 Accessories

Accessory	Order No.
O ₂ sensor master InPro 6950	52 206 113
O ₂ sensor simulator	52 200 891
Cap sleeve with protective cage (P type)*	52 200 038
Cap sleeve w/o protective cage (N type)*	52 200 037
VP cable VP6-ST/3 m	52 300 108
VP cable VP6-ST/5 m	52 300 109
VP cable VP6-HT/3 m	52 300 112
VP cable VP6-HT/5 m	52 300 113
VP cable VP8-ST/3 m	52 300 354
VP cable VP8-ST/5 m	52 300 355
VP cable VP8-HT/3 m	52 300 361
VP cable VP8-HT/5 m	52 300 362
Oxygen zeroing gel (3 × 25 ml)	30 300 435

* 1.4404

For other cable lengths or types, please contact your local METTLER TOLEDO representative.

9.3 Spare parts

Spare part	Order No.
Membrane body, single, InPro 6950	52 206 105
Membrane kit InPro 6950 (4 membrane bodies, 1 O-ring set, 2 × 5 ml electrolyte bottle, watering cap)	52 206 106
Membrane body, single, InPro 6950 Gas	52 206 108
Membrane kit InPro 6950 Gas (4 membrane bodies, 1 O-ring set, 2 × 5 ml electrolyte bottle, watering cap)	52 206 109
InPro 6950 electrolyte pack (3 × 5 ml)	30 298 426
Interior body InPro 6950 (with quick disconnect)	52 206 112

9.4 Recommended transmitters

Transmitter	Order No.
O ₂ transmitter M 700, coated C	52 121 171
O ₂ transmitter M 700, stainless steel S	52 121 174
O ₂ transmitter M 700, coated, Ex VPW	52 121 172
O ₂ transmitter M 700, coated, Ex 24V	52 121 173
O ₂ transmitter M 700, stainl. steel, Ex VPW	52 121 175
O ₂ transmitter M 700, stainl. steel, Ex 24V	52 121 176

9.5 Recommended housings

Housing (12 mm Ø)	Order No.
Static housing	
InFit 761 CIP	
Retractable housings	Please ask your local METTLER TOLEDO representative.
InTrac 777e	
InTrac 797e	
Immersion housing	
InDip™ 550	

Note: The housings are available in different versions.

Please contact your local METTLER TOLEDO sales

organization to get the right ordering information.

10 Theory of the polarographic sensor

10.1 Introduction

Two types of electrodes are employed in analytical work: **potentiometric** and **amperometric** electrodes.

- Potentiometric electrodes develop a voltage generated by the activity of a particular ion. Examples of such electrodes are glass electrodes (like pH electrodes) and most ion-selective electrodes. Their individual potentials cannot be determined. The measurable quantity is the difference of potential between the measuring electrode and an inert reference electrode. The potential of the reference electrode must be constant.

All potentiometric electrodes are subject to Nernst's law and for this reason electrodes and measuring instruments are in most cases interchangeable.

An important requirement of potentiometric measurements is the virtually currentless determination of the electrode voltage. During measurement no chemical reaction occurs and the solution remains in equilibrium.

- In the case of **amperometric electrodes**, such as the **oxygen electrode**, activity measurement is based on a current measurement.

The conventional oxygen electrode consists of a cathode and an anode conductively connected by an electrolyte. A suitable polarization voltage between the anode and the cathode selectively reduces the oxygen at the cathode.

Cathode reaction:

Anode reaction:

In the case of the InPro 6950 sensor, the measurement system has been expanded by reference and guard ring electrodes. This design modification allows the precise concentration of oxygen to be measured down to trace level. In contrast to conventional sensors, this also changes the electrochemical reaction at the anode..

InPro 6950 anode reaction:

These chemical reactions result in an electric current which is proportional to the oxygen partial pressure (pO₂). The oxygen electrode consumes oxygen which is continuously extracted from the solution. The viscosity and flow rate of the solution are therefore important parameters.

The electrode current of an oxygen electrode is determined not only by the oxygen partial pressure but by many other electrode parameters. The electrode currents of different electrode types may differ by several powers of ten. For this reason oxygen electrodes and amplifiers cannot be freely interchanged.


10.2 Principle of the design of an oxygen electrode

There are two main types of oxygen electrodes:

- Electrodes **without** membrane
- Electrodes **with** gas-permeable membrane (Clark Principle)

The membrane electrode according to Clark is today most widely used. As compared to the electrode without membrane it possesses the following advantages:

- Oxygen measurement in gases and solutions
- No mutual contamination of electrode and solution
- No or little dependency on flow

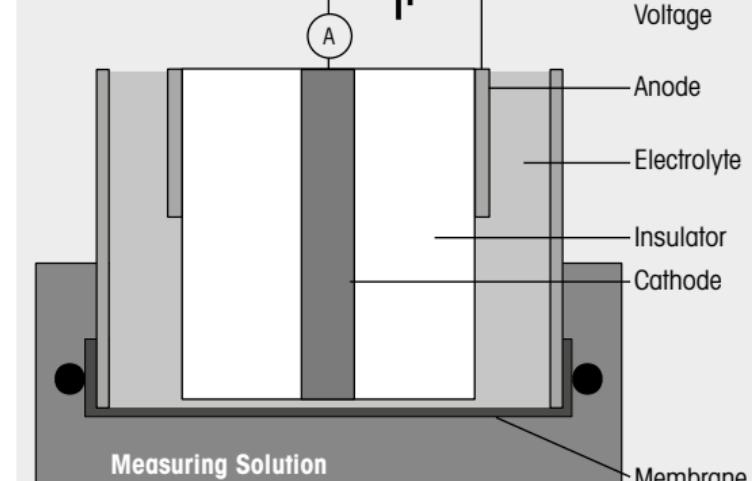
In the case of the Clark electrode, geometrical configuration is very important. In particular, the thickness of the electrolyte film between the cathode and the membrane must be within narrow tolerances so as to ensure good linearity and a low zero current (current in a nitrogen atmosphere). The figure above exemplifies the principle of the design of a Clark-type oxygen electrode.

Principle of the design of the InPro 6950

The InPro 6950 sensor is based on the Clark electrode. However it also has a cathode ring (guard ring electrode) and a reference electrode.

The guard ring electrode is connected to a separate polarization voltage and hence consumes the oxygen that is released from the electrolytes to the measurement cathode.

The reference electrode is used as a reference point to adjust the necessary polarization voltage between the cathode and anode and hence creates stable conditions within the sensor during measurement. This guarantees exceptional signal stability in the measurement system down to trace levels of oxygen.

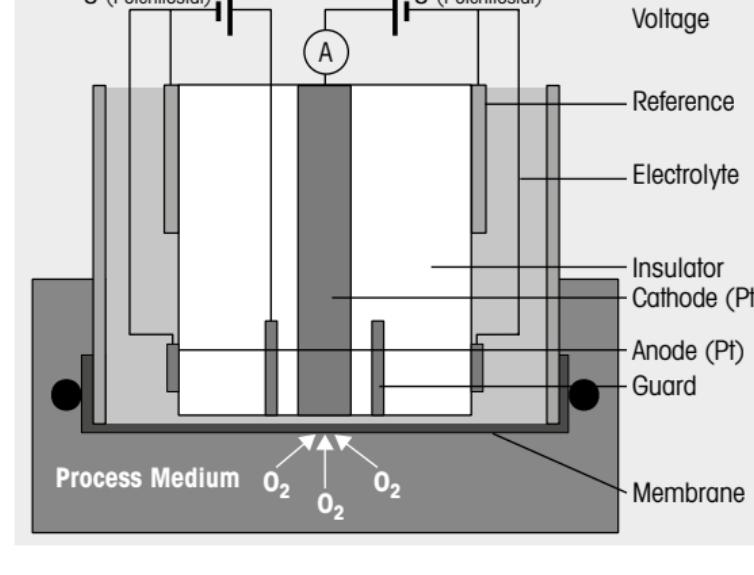10.3 Parameters determining current

The quantity of oxygen diffused in and the magnitude of the electrode current are influenced by the following parameters:

- Oxygen partial pressure of the solution
- Membrane material and thickness
- Size of cathode
- Polarization voltage
- Temperature
- Flow conditions in the solution

Fick's law gives the mathematical relationship between these parameters:

$$I = k \times D \times a \times A \times \frac{pO_2}{x}$$

I = Electrode current

k = Constant

D = Coefficient of diffusion of O₂ in the membrane

a = Solubility of O₂ in the membrane material


A = Cathode surface

pO₂ = Oxygen partial pressure in the solution

x = Thickness of gas-permeable membrane

10.4 Polarization voltage

The voltage between the anode and the cathode is so selected that the oxygen is fully (> A, see polarogram) reduced while other gases are unaffected (< D). The ideal voltage for the Pt/Ag/AgCl or Pt/Pt/Ag/AgCl system is between -500 and -750 mV.

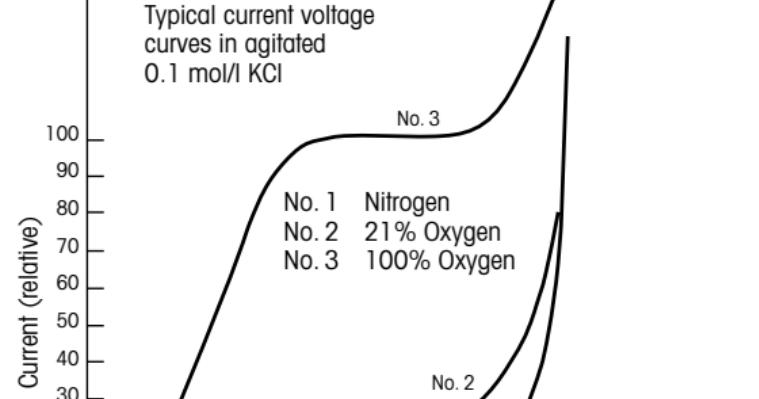The polarization voltage should remain as constant as possible. Besides a constant source of voltage the following prerequisites must be met: The electrical resistance of electrolyte film must not exceed a particular value so as to avoid a voltage drop.

The anode must have a large surface so as to prevent polarization of the anode by the electrode current.

10.5 Temperature

The temperature dependence of the current passing through an oxygen electrode when referred to a constant O₂ partial pressure is determined mainly by the properties of the gas-permeable membrane.

10.6 Dependence on flow

With most oxygen electrodes the electrode current is smaller in stagnant solutions than in agitated ones. The consumption of oxygen by the electrode results in extraction of oxygen from the solution in the close vicinity of the cathode outside the membrane. The oxygen is replaced by diffusion. If the electrode current is strong, the solution cannot fully restore the oxygen by diffusion. This results in an electrode current weaker than would correspond to conditions in the solution. In agitated solutions the oxygen is conveyed to the surface of the membrane not only by diffusion but additionally by the flow (convection). In that case no oxygen impoverishment occurs at the membrane surface.

A high degree of flow dependence occurs mainly with large cathodes, thin and highly permeable membranes, i.e. where electrode currents are large.

The problem of flow dependence is often solved by prescribing a minimum flow rate.

In METTLER TOLEDO InPro 6950 electrodes, the thin PTFE membrane determining the electrode current (i.e. the actual measuring signal) is separated from the sample solution by a relatively thick silicone membrane. This latter is highly permeable to oxygen molecules and thus acts as an oxygen reservoir. The diffusion of oxygen out of the sample solution into the silicone membrane is spread over a wide area. Since this results in less oxygen being extracted from the sample solution per unit area, the double PTFE/silicone membrane forms an effective buffer against disturbances due to hydrodynamic flow.

This membrane coupled with the guard ring and the special electrolyte ensures excellent signal stability even when the hydrodynamic flow stops (as on a beer filling line for example).

10.7 Oxygen partial pressure – oxygen concentration

The electrode current is dependent on the partial pressure of oxygen nad the oxygen permeability of the membrane. Conversion of partial pressure into concentration of oxygen depends on the measurement medium (measurement in liquids or gases).

Measurement in liquids

When measuring in liquids, the concentration of oxygen also depends on the solubility of the oxygen in the measurement medium. However, since this is not measured by the sensor current, the concentration of oxygen must be calculated in the transmitter. To do this, Henry's law is applied which states that the concentration of oxygen is proportional to the partial pressure of oxygen (pO_2).

$$Cl = pO_2 \times a$$

a = Solubility factor

If "a" is constant, the oxygen concentration can be determined by means of the electrode. This applies at constant temperature and with dilute aqueous solutions such as drinking-water.

The solubility factor is strongly influenced not only by the temperature but also by the composition of the solution:

Medium, saturated with air	Solubility at 20 °C (68 °F) and 760 mm Hg
Water	9.2 mg O ₂ /l
4 mol/l KCl	2 mg O ₂ /l
50 % Methanol-water	21.9 mg O ₂ /l

Although the solubilities vary widely, the oxygen electrode gives the same reading in all three solutions.

Thus, determination of the oxygen concentration is only possible with constant and known solubility factors "a".

Solubility may be determined by a Winkler titration or the method developed by Käppeli and Fiechter.

Measurement in gases

The concentration of oxygen when measuring in gases is always given as a proportion by volume of the gas composition. Common units are % (by vol.) ppm (by vol.).

They can simply be converted from one unit of measurement to the other.

Example:

The percent by volume of the composition of air is generally known. For example, air contains 20.95 % oxygen. This corresponds to 209,500 ppm (by vol.). (Conversion: ppm value = 10,000 × value in %)

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US-Normen

Titel	Nummer	Ausgabe-Datum
Approval Standard for Electrical Equipment for Use in Hazardous (Classified) Locations – General Requirements	FM Class 3600	2022
Approval Standard for Intrinsically Safe Apparatus and Associated Apparatus for Use in Class I, II & III, Division 1, Hazardous (Classified) Locations	FM Class 3610	2021
Approval Standard for Electrical Equipment for Measurement, Control and Laboratory Use	FM Class 3810	2021
Explosive Atmospheres – Part 0: Equipment – General Requirements	ANSI/UL 60079-0:2020	2020
Explosive Atmospheres – Part 11: Equipment Protection by Intrinsic Safety "i"	ANSI/UL 60079-11:2018	2018
Safety Requirements for Electrical Equipment for Measurement, Control, and Laboratory Use – Part 1: General Requirements	ANSI/UL 61010-1:2018	2018

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Titel	Nummer	Ausgabe-Datum
Explosive Atmospheres – Part 0: Equipment – General Requirements	CAN/CSA-C22.2 No. 60079-0:2019	2019
Explosive Atmospheres – Part 11: Equipment Protection by Intrinsic Safety "i"	CAN/CSA-C22.2 No. 60079-11-14:R2018	R2018
Safety Requirements for Electrical Equipment for Measurement, Control, and Laboratory Use – Part 1: General Requirements	CAN/CSA-C22.2 No. 61010-1-12:R2017	R2017

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Titre	Numéro	Date de publication
Approval Standard for Electrical Equipment for Use in Hazardous (Classified) Locations – General Requirements	FM Class 3600	2022
Approval Standard for Intrinsically Safe Apparatus and Associated Apparatus for Use in Class I, II & III, Division 1, Hazardous (Classified) Locations	FM Class 3610	2021
Approval Standard for Electrical Equipment for Measurement, Control and Laboratory Use	FM Class 3810	2021
Explosive Atmospheres – Part 0: Equipment – General Requirements	ANSI/UL 60079-0:2020	2020
Explosive Atmospheres – Part 11: Equipment Protection by Intrinsic Safety "i"	ANSI/UL 60079-11:2018	2018
Safety Requirements for Electrical Equipment for Measurement, Control, and Laboratory Use – Part 1: General Requirements	ANSI/UL 61010-1:2018	2018

Normes canadiennes

Titre	Numéro	Date de publication
Explosive Atmospheres – Part 0: Equipment – General Requirements	CAN/CSA-C22.2 No. 60079-0:2019	2019
Explosive Atmospheres – Part 11: Equipment Protection by Intrinsic Safety "i"	CAN/CSA-C22.2 No. 60079-11-14:R2018	R2018
Safety Requirements for Electrical Equipment for Measurement, Control, and Laboratory Use – Part 1: General Requirements	CAN/CSA-C22.2 No. 61010-1-12:R2017	R2017

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All rights reserved